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
Artikelnummer: 653 375

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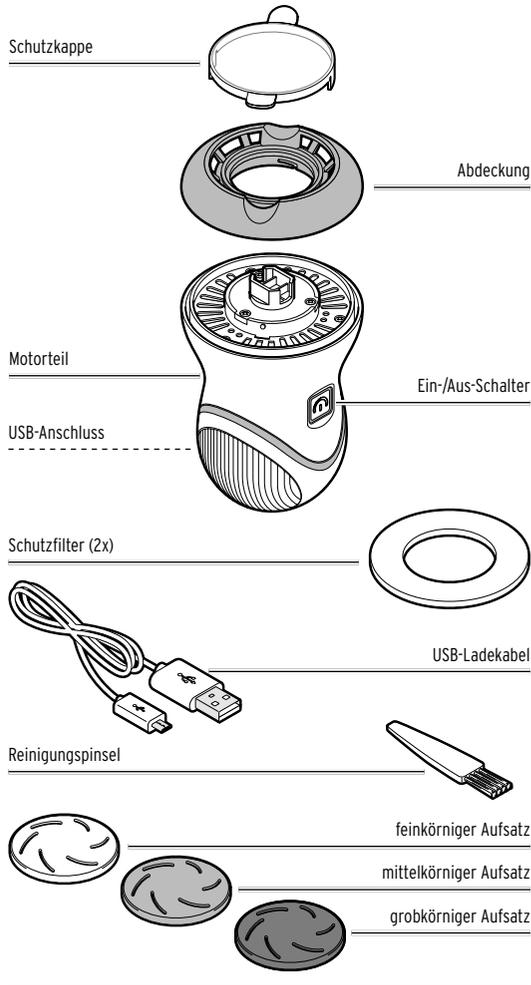
At a glance (contents)

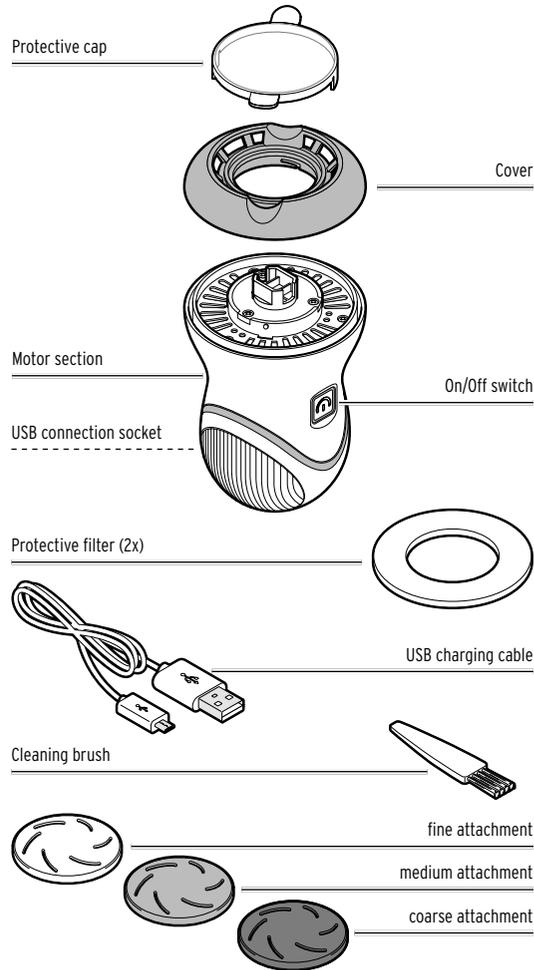Safety warnings

This product is equipped with safety features. Nevertheless, read the safety warnings carefully and only use the product as described in these instructions to avoid accidental injury or damage.

Keep these instructions for future reference. If you give this product to another person, remember to also include these instructions.

Intended use

- The device is intended for removing hard skin on feet, hands and elbows.
- The product is designed for domestic use and is not suitable for commercial purposes, and especially not for medical foot care.

Do not use the device...

... if you suffer from diabetes or other wound-healing disorders. The possibility of injuring yourself with this device is very slight. Nevertheless, we are obliged to draw your attention to the fact that you should not use it if you suffer from the illnesses mentioned above.

... on reddened or damaged skin or on burns. Other areas without hard skin, as well as birthmarks, moles, warts and other skin features must not be treated with the product.

... if another person has already used it. The device and attachments that come with it should be used by one person only. This will prevent the spread of infections such as athlete's foot.

Danger to children and to people with limited ability to operate devices

- Children must not be allowed to play with the device. Keep children under 8 years of age away from the device.
- Keep the packaging materials out of the reach of children. They pose several risks, including the risk of suffocation!
- This device must not be used on infants, as their skin is very sensitive.
- This device may be used by children aged 8 and over and by people with diminished physical, sensory or mental capacity or lacking the requisite knowledge/experience, provided they are supervised or have been instructed in the safe use of the device and understand the potential dangers.
- Children must not clean the device unless they are over 8 years of age and are supervised when doing so. The device is maintenance-free.

Danger due to explosion

- The device contains a built-in lithium battery. Do not throw the device into a fire. Risk of explosion!

Risk of injury

- Stop removing hard skin immediately if you start to feel pain!
- Do not hold the device/attachment on one area, as the friction causes intense heat to develop.
- Keep long hair away from the rotating attachment during use.

Material damage

- Protect the device from moisture. Do not immerse it in water or other liquids.
- The rechargeable battery is firmly built into the device and you must not attempt to replace it by yourself. If damage occurs to the rechargeable battery, consult a specialist workshop.
- Protect the device from hard impacts, falls, dust, moisture, direct sunlight and extreme temperatures.
- Stop using an attachment if it is very worn down. Do not use the device if it is damaged. Only have repairs to the device carried out by a specialist workshop.
- Clean the device thoroughly after each use as described in the chapter "Cleaning and disinfecting the device" to remove any grease deposits and other residues.
- Do not use any caustic, aggressive or abrasive products for cleaning.
- When operating or charging the rechargeable battery, ensure that ... you observe the technical specifications of the mains adapter before connecting the device to the power supply using a USB mains adapter. The device must be connected to a mains adapter with protection class II (symbol: .


Please note

Hard skin is normal, especially on feet. However, it should not become too thick or hard as it could then become painful. Excessive formation of hard skin can be caused by shoes that are either too small or have an improper fit, incorrect foot positioning or extremely dry skin. If you are suffering from excessive hard skin, try to fix the causes and consult a doctor for advice.

To prevent the excessive formation of hard skin, regular foot baths are recommended, followed by rubbing your feet with lotion or foot creams. Make sure that your shoes fit well and try to wear different shoes every day if possible.

Charging the battery

The rechargeable battery must be fully charged prior to first use. The device must be switched off to do so. **Do not use the device while it is being charged.**

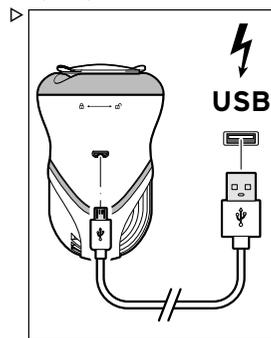Connect the device to a mains adapter with a USB connection socket (not included), as illustrated. The device will start to charge automatically. The indicator light will light up in red. Charge the rechargeable battery until the indicator light turns green. Charging time for a flat battery: approx. 150 minutes
Operating time when the battery is fully charged:
approx. 40 minutes (speed setting 1),
approx. 38 minutes (speed setting 2).

Use

CAUTION - risk of injury

- Stop removing hard skin immediately if you start to feel pain!
- Keep long hair away from the rotating attachment during use.

Always leave a little bit of hard skin remaining so that the skin is protected.

Changing attachments


-  The coarse attachment is intended for removing thicker layers of hard skin.
-  The medium attachment is intended for touching up remaining hard skin after the coarse attachment has been used.
-  The fine attachment is intended for removing thinner layers of hard skin.

1. Switch off the device.
2. Check the device for signs of damage before every use. Do not use any damaged or severely worn attachments.

To remove the attachment, use the spaces on either side to grasp the attachment and pull it up.


To connect the attachment, place it onto the device as illustrated and push it down as far as it will go.

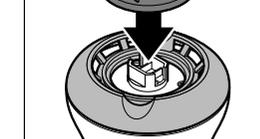
Inserting the protective filters

Unscrew the cover anticlockwise towards  and lift it up to remove it.

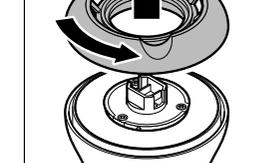Place a filter in the slot and screw the cover back on towards .

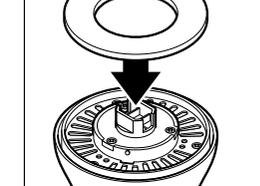Removing hard skin


CAUTION - risk of injury

Do not hold the attachment on one area, as the friction causes intense heat to develop.

We recommend soaking your feet for approx. 10 minutes in order to soften your skin.

1. Use your hand to stretch the area of the skin you wish to treat.

2. **To switch the device on,** hold down the On/Off switch  for approx. 3 seconds. The indicator light in the On/Off switch lights up green. The attachment begins to rotate at a normal speed (speed setting 1).

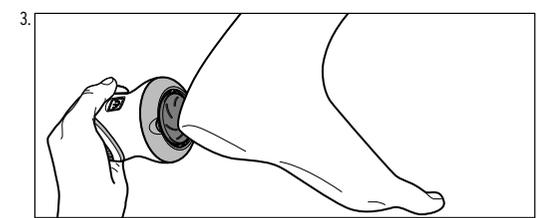While applying gentle pressure, carefully move the rotating attachment over the hard skin you wish to remove.

- ▷ **To increase the speed,** briefly press the On/Off switch once more. The indicator light in the On/Off switch lights up red. The attachment begins to rotate at a higher speed (speed setting 2).

4. **To switch the device off,** hold down the On/Off switch  for approx. 3 seconds. The indicator light on the underside goes out.

The device is equipped with a protection mechanism and will automatically switch itself off after 5 minutes. If you want to continue the treatment, simply switch the device on again.

After treatment

- ▷ Rub a suitable body cream into the area where the skin was treated. Do not use any care products that might irritate the skin, e.g. deodorant.
- ▷ After use, clean the attachments as described in the section titled "Cleaning and disinfecting the device".
- ▷ Always place the protective cap onto the attachment after use.

Cleaning and disinfecting the device

- ▷ Clean the attachments using the cleaning brush provided.
- ▷ If necessary, wipe the device with a cloth that has been slightly dampened in water.
- ▷ Disinfect the attachments, as required, using ethyl alcohol available to buy at a chemist.

Technical specifications

Model:	653 375
Input:	5 V / 2 A micro-USB
Rechargeable battery	Lithium-ion 3.7 V DC 450 mAh, Watt-hour rating: 1.67 Wh (tested as per UN 38.3) (non-replaceable)
Protection class:	III 
Ambient temperature:	+10 to +40 °C
Manufactured by:	Commaxx B. V. Wiebachstraat 37 6466 NG Kerkrade The Netherlands

In the course of product improvement, we reserve the right to make technical and optical modifications to the product.

Disposal

This **product**, its **packaging** and the built-in **rechargeable battery** have been manufactured from valuable materials that can be recycled. Recycling reduces the amount of refuse and helps to preserve the environment.

Dispose of the packaging at a recycling point which sorts materials by type. Make use of the local facilities provided for collecting paper, cardboard and lightweight packaging.

 **Devices** marked with this symbol must not be disposed of along with normal household refuse! You are legally obliged to dispose of old devices separately from household waste. Electronic devices contain hazardous substances and, if stored or disposed of improperly, may cause harm to health and the environment. Information about collection points where old devices can be disposed of free of charge is available from your local authorities.

 Single-use batteries and rechargeable batteries must not be disposed of along with household waste! You are legally obliged to dispose of flat single-use batteries and flat rechargeable batteries at a local authority collection point or to return them to a battery retailer.

 Warning! This device contains a battery that, for safety reasons, is firmly built into the device. It cannot be removed without destroying the casing. Improper disassembly leads to safety hazards. For this reason, take the device intact to a collection point where the device and battery will be disposed of properly.